

Zygoma Kit is the first ultrasonic tips' kit for zygomatic implants

ADVANTAGES

- > 4 calibrated tips , **non sterile**
- > more agile, effective and safe surgical approach
- > safe and predictable protocols: in addition to the absolute precision of positioning the ultrasonic Tips soft tissues are preserved .

FEAUTURES AND SIZES

Inserts' diameters:

Zigoma Kit	Ø	L
ES052XZT	3,0 mm	74 mm
ES052XZLT	3,5 mm	94 mm
ES3.5ZT	3,0 mm	74 mm
ES3.5ZLT	3,5 mm	94 mm

- > Marking with depth's notches :
- > Long inserts: 30-35-40-45-50-55 mm
- > Short inserts: 30-35-40-45-mm

WARNING:

High Efficiency's tips: it's mandatory only one use, to avoid the risk of breakage.

Esacrom declines any responsibility in case of reuse.

Only to be used with ESACROM'S medical devices

For any further information please contact:

Esacrom srl
VIA ZAMBRINI, 6/A – 40026 IMOLA(BO)
TEL. +39 0542 643527
FAX +39 0542 482007
Web-site: www.esacrom.com



Dr. Marco Mozzati

Bachelor of Medicine and Surgery in 1986 and specialized in Dentistry at the University of Turin in 1990., in 1998 he was Medical Director at the first level by the Autonomous Service of Oral Rehabilitation and Prosthetics Maxillofacial and Dental Implants of the St. John the Baptist (TO) in 2002: dental consultant at the hospital CTO, Department of Plastic Surgery Traumatology and Maxillofacial Prof. Magliacani. In 2005 Medical Director of the Head SSCVD Stomatology Surgery Department of Dentistry directed by Prof. Stefano Carossa, c / o the Hospital San Giovanni Battista (TO). Adjunct professor since 2004 with the teaching of Oral Surgery for the Degree Course in Dentistry. Since 2005 the teaching of Oral Surgery for the Degree Course in Dental Hygiene at the University of Turin. Founding Member of Piezosurgery Academy. Active member SICOI. Founding member and President of ANTHEC (Academy of transfusional not Hemo-Components). Author of scientific articles and publications regarding oral surgery and implantology.



Dr. Roberto Pistilli

Bachelor of Medicine and Surgery of the University of Rome in 1983 Specialist in Maxillofacial Surgery since 1988; executive Maxillofacial surgeon at the ACO San Camillo in Rome from 1989 to 1994, at the ACO S. Giovanni Addolorata in Rome from 1994 to 1999 and at the ACO San Filippo Blacks (Rome) from 1999 to present with appointment of highly professional (P1). Holder of various courses of Oral and Maxillofacial National and International. 2007-2008: Assistant Professor in the years at the Faculty of Dentistry at the University of Guarulhos (São Paulo), Brazil, and in 2007 at the Faculty of Medicine and Dentistry of New Jersey University (UMDNJ). Director of the Annual Course of Specialization in Oral Surgery and Implantology at ACO San Filippo Blacks (Rome) - Professor in the Course on "Microscopy operating in oral surgery" - The University G. D'Annunzio - (CH) Dir. Prof. A. Piattelli. 2010/2011: Adjunct Professor in the Master of Implantoprosesi at La Sapienza University (Rome). 2011/2012 at the University Alma Mater Studiorum - (BO) Dir. Prof. L. Checchi. Active Member SIO.- Member of the Scientific Committee of the international journal "Journal of Osteointegration", Editorial Board of the journal "European Journal of Oral Implantology." President of A.R.C.O.I. Author of over 70 publications in national and international scientific and speaker at numerous courses and conferences in Italy and abroad. Collaborator for several scientific expansions (crestal). He is interested mainly in regenerative surgery of the jaw and advanced implantology.



Dr. Pietro Felice

Dr. Felice graduated in Dentistry and Dental Implants and later in Medicine and Surgery at the Alma Mater Studiorum, University of Bologna. He later completed his PhD in Dermatological Sciences and Reconstructive Maxillofacial and Facial Plastic. He is currently a researcher at the Alma Mater Studiorum of Bologna. In the same University he teaches in Implant Surgery Degree in Dentistry. He is the author of numerous scientific articles in international journals. He is a member of the editorial board of international journals and a member of the Cochrane Collaboration. Exercising professional activity limited to oral surgery and implant.

ZYGOMA KIT

ULTRASONIC DRILL

*THE ONE and ONLY
ULTRASONIC INSERTS KIT FOR
ZYGOMATIC IMPLANTS*






Surgysonic®

18/05/16



Clinical case

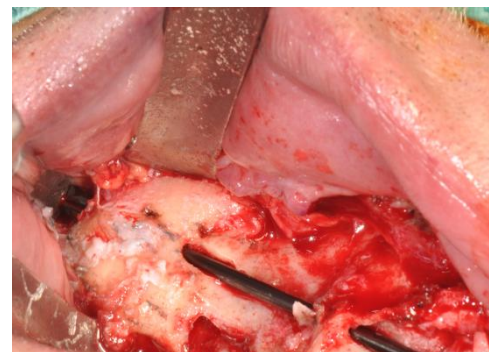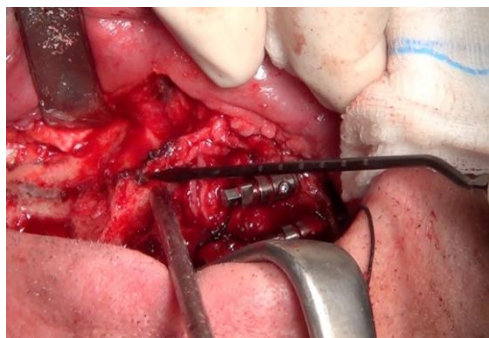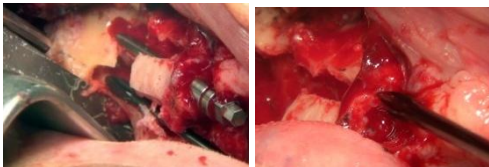

**INSERTI PER
ULTRASONIC DRILL**



ES052XZT	ES052XZLT	ES3.5ZT	ES3.5ZLT
L. 74 mm	L. 94 mm	L. 74 mm	L. 94 mm
Ø 3,0 mm	Ø 3,5 mm	Ø 3,0 mm	Ø 3,5 mm

Recommended parameters for inserts ' proper use with Surgysonic®

Power	(U)	50
Vibration	(V)	80
Max Power *(U)		60

•Maximum power for proper insert use

•**WARNING:** With high power high risk of breakage. **Mandatory: only 1 use**

Clinical case

